

INTELLIGENT[®] SHAVER 2500



INTELLIGENT[®] SHAVER 2500

Electric Power Shaver System for removal
of soft and hard tissue during Arthroscopic procedures.

Intelligent Shaver 2500 with **Opti-Cut[™]** Algorithm for optimum and superior resection of soft tissue. Faster, simpler, trendier with the maximum torque in its class equipped with the latest opti-cut [™] Algorithm enabling more power for effective control and more efficient bone debriement and tissue resection.

SPECIAL FEATURES AND BENEFITS

- Utilizing a unique Computerised Servo Controller Technology, resulting in unmatched power.
- Complete autoclavable system.
- Improved control and flow of suction line.
- Foot switch for Right, Left and Oscillation control.
- Operating speed upto 5500 rpm.
- Integrated self test and diagnostic read out.
- MDE incorporates the state-of-art high Performance 3 Ø BLDC motors.
- Standard high speed shaver hand piece with rotation speed of 600 rpm to 5500 rpm.
- Ergonomically designed foot switch with dedicated Soft touch pedals that stays clean & is fully immersible.
- Wide screen graphic display of all operating parameters.
- Oscillation speed from 600 rpm up to 3500 rpm.
- Ergonomically designed lightweight handpiece with easy to use switch for suction control.

MDE




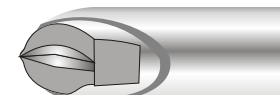


INTELLIGENT[®] SHAVER 2500

TECHNICAL SPECIFICATIONS

Supply	: 230 V Ac, 50 Hz.
Control	: Micro computerized servo controller
Speed Ranges	: In 6 stages of 600, 1000, 2200, 3500, 4500, 5500
Frequency	: 50 / 60 Hz.
Electro servo drive	: Hand piece with constant torque BLDC motor.
Sterilization	: Hand piece - Soakable / auto clavable
Foot switch	: Sokable / water proof
Classification according to MDD	: CE
Self Diagnosis	: Self diagnosis check for all important functions and all accessories.

Note : Specifications are subject to change without notice

ORDERING INFORMATION

Arthroscopy Shaver system		LS 100 00000
Accessories		
Shaver Handpiece Electro Servo Drive		LS 100 20000
Foot Switch		LS 100 30000
Cable 3 mtrs.		LS 100 40000
Blades		
Full Radius Resector		LS 100 50001
Aggressive Full Radius		LS 100 50002
Meniscus Cutter		LS 100 50003
Round Burr		LS 100 50004
Tapered Burr		LS 100 50005
Oval Burr		LS 100 50006